

Advanced Venepuncture Arm User Guide



Advanced Venepuncture Arm - Light Part No: **00290**
Advanced Venepuncture Arm - Brown Part No: **00298**
Advanced Venepuncture Arm - Black Part No: **00296**



Designed and manufactured by Limbs & Things Limited,
Sussex Street, St. Philips, Bristol, BS2 0RA, UK.
Telephone: +44 (0)117 311 0500 Fax: +44 (0)117 311 0501
sales@limbsandthings.com www.limbsandthings.com

Issue 1, July 2008
© 2008 Limbs & Things

Contents

Product information

Components

READ THIS FIRST: Basic fluid unit operations 1

Getting the Arm ready to use

Filling the fluid unit 3

Priming the Veins 4

Using the Arm 5

Emptying the Arm

Emptying the veins prior to replacement 7

Emptying the fluid unit prior to storage/transport 8

Removing the Skin 9

Fitting the Skin 10

Changing the Veins 11

Storing the Arm 13

Advanced Venepuncture Arm - Light

Part No: **00290**

This product provides the ideal platform for all healthcare trainees needing to practise the skills related to venepuncture and IV cannulation. Designed with durability, realism and ease of use in mind, this product addresses the fundamental needs of clinical skills tutors, trainees and technicians alike.

Skills

- Venepuncture
- IV Cannulation

Features

- Veins are:
 - leak resistant when used at the correct pressure
 - easy to replace
- Can be used with vacuum blood collection systems
- Realistic blood flashback from pressurised fluid system
- Realistic, soft, flexible skin and underlying palpable veins
- Easy to clean, service and maintain
- Accessible veins: dorsal metacarpal, cephalic, basilic, and median cubital are easy to replace
- 'Closed' blood flow system with no messy bottles to change
- Pressure of blood can be easily adjusted
- Can be adapted for professional-to-patient communication (role play) using **AV Arm Harness** Part No: **00300**

Package supplied

- 1 **AV Arm Skin - Light** Part No: **00291**
- 1 **AV Arm Veins** Part No: **00292**
- 1 **AV Vein Module** Part No: **00295**
- 1 **AV Mock Blood (500ml)** Part No: **00294**
- 1 **Refill Bottle** Part No: **00303**
- 1 **AV Arm Rest** Part No: **00302**
- 2 lubricant sachets
- 1 arm

Advanced Venepuncture Arm - Brown

Part No: **00298**

As per the Advanced Venepuncture Arm - Light.

Package supplied

- 1 **AV Arm Skin - Brown** Part No: **00299**
- 1 **AV Arm Veins** Part No: **00292**
- 1 **AV Vein Module** Part No: **00295**
- 1 **AV Mock Blood (500ml)** Part No: **00294**
- 1 **Refill Bottle** Part No: **00303**
- 1 **AV Arm Rest** Part No: **00302**
- 2 lubricant sachets
- 1 arm

Advanced Venepuncture Arm - Black

Part No: **00296**

As per the Advanced Venepuncture Arm - Light.

Package supplied

- 1 **AV Arm Skin - Black** Part No: **00297**
- 1 **AV Arm Veins** Part No: **00292**
- 1 **AV Vein Module** Part No: **00295**
- 1 **AV Mock Blood (500ml)** Part No: **00294**
- 1 **Refill Bottle** Part No: **00303**
- 1 **AV Arm Rest** Part No: **00302**
- 2 lubricant sachets
- 1 arm

Important: handling & hygiene

- Venepuncture training is a potentially messy process so always set the product up on a wipe clean surface. The instructions supplied in this User Guide are designed to minimise the risk of any spillage. Even so, it is a good idea to keep paper towelling close to hand just in case.
- Do not wear 'good' clothes when training with this product. If mock blood does get onto clothing it will not permanently stain them and can be washed out using a biological detergent.
- If mock blood gets into the eyes flush with plenty of cold water.
- Observe normal hygiene procedures after handling the product and mock blood.

Components



AV Arm Skin - Light

Part No: **00291**

(Depending on version bought)



AV Arm Skin - Brown

Part No: **00299**

(Depending on version bought)



AV Arm Skin - Black

Part No: **00297**

(Depending on version bought)



AV Arm Rest

Part No: **00302**



AV Arm Veins

Part No: **00292**

(With adhesive fastening strips)



AV Mock Blood (500ml)

Part No: **00294**

(Replacement part)



Water-based Lubricant

Part No: **00293**



Fluid Unit

(Not user serviceable. Must be returned to Limbs & Things for service or repair)





Refill Bottle

Part No: **00303**


READ THIS FIRST: Basic fluid unit operations

Connecting

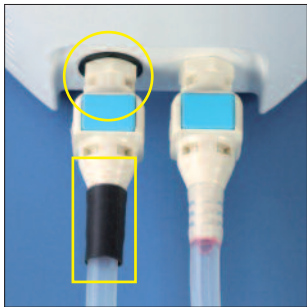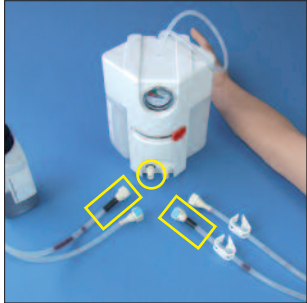

When connecting the refill bottle or the Arm to the fluid unit, ensure that the tube with the black sleeve attaches to the left hand connector with the black ring on it.

Selecting

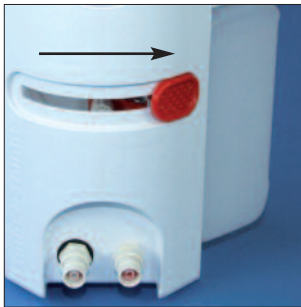


The red reservoir selector lever on the front of the fluid unit selects which of the reservoirs will be pressurised.

Slide it to the left if you intend to pressurise the left reservoir.

Slide it to the right if you intend to pressurise the right reservoir.


Pressurising

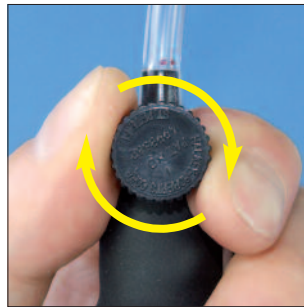

Turn the release valve clockwise to close it. Then pressurise the fluid unit by squeezing the pressure bulb.

NOTE
There is a pressure release valve inside the fluid unit. This prevents damage caused by overpressurisation. If the cutoff pressure is reached a hissing sound will be heard as the air is vented and the pressure is dropped to a safe level.


Depressurising

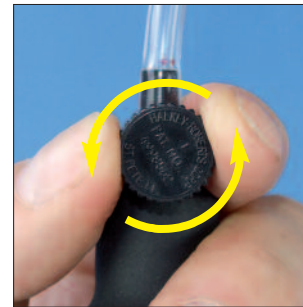

Turn the release valve anti-clockwise to open it when you need to depressurise the fluid unit.

Filling the fluid unit

1

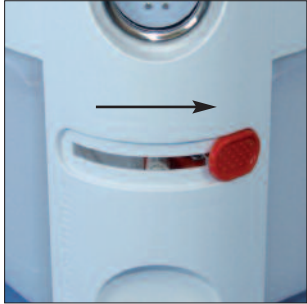


IMPORTANT
When filling the fluid unit, blood must always be fed into the left hand reservoir.

Slide the reservoir selector to the right.

Connect the refill bottle.

2

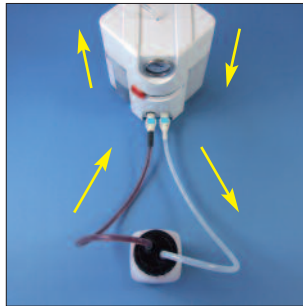

Close the valve on the pressure bulb. Pressurise the fluid unit into the grey zone.

Air in the right reservoir flows into the bottle forcing the blood out of it and into the left hand reservoir.

Keep the pressure in the grey zone to ensure that the reservoir fills quickly.

3



Once the blood reaches the maximum level of 300ml in the left reservoir, open the pressure valve to depressurise the fluid unit and to stop blood transferring.



Disconnect the refill bottle.

The fluid unit is full and ready to be attached to the Arm.

3

Priming the Veins

1

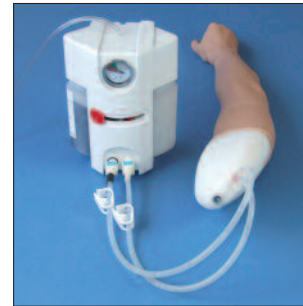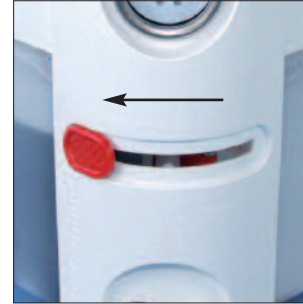

Having carried out 'Filling the fluid unit' (page 3), slide the reservoir selector to the left.

Connect the Arm.

2




Ensure both of the white clips on the arm tubes are open.

Close the valve on the pressure bulb. Pressurise the fluid unit into the green zone.

IMPORTANT
Do not overpressurise (into the grey zone) as this will cause blood to leak from the puncture sites under the Skin.

3

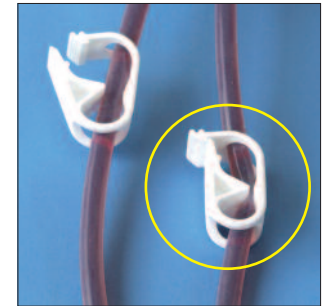

Air pressure in the left reservoir forces the blood out through the Arm and into the right reservoir.

When there are no more air bubbles coming out of the Arm close the white clip on the right hand tube. This will maintain the correct pressure in the Arm during training.

The Arm is ready to use.

4

Using the Arm

Taking blood



Blood can be taken from the arm by using vacuum collection devices such as Vacutainer, or the more traditional needle and syringe method.

IV cannulation can also be performed on the product.


Maintaining pressure



As the blood level decreases in the left reservoir the pressure will drop.



For best results maintain the pressure in the green zone during training.

Removing trapped air

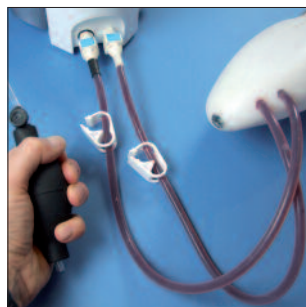

Occasionally release the clip on the right hand tube. This will allow blood to flow through the Arm to remove any air bubbles that might have been introduced during training.


Reversing flow

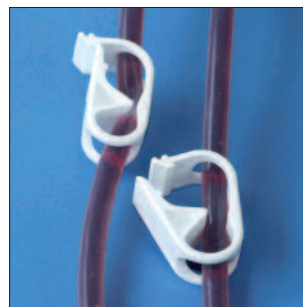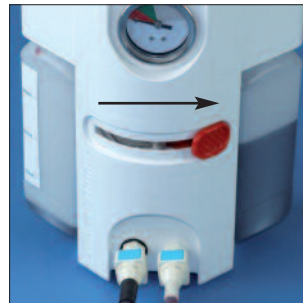

The Arm can also be pressurised from the right reservoir when there is a sufficient volume of blood in it.

Slide the reservoir selector to the right. Open the clip on the right tube and close the clip on the left tube. Ensure the pressure is in the green zone.

Refilling the bottle


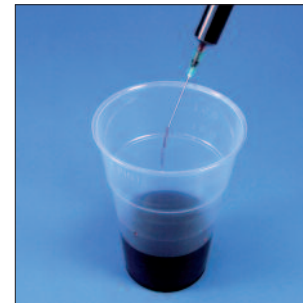


To reuse the blood squirt the contents of the syringe into a container. At the end of the training session the blood can be poured back into the refill bottle.

Follow **'Filling the fluid unit' (page 3)** when the left reservoir is empty.

Emptying the veins prior to replacement

1

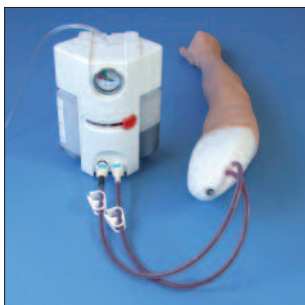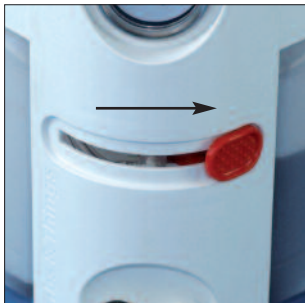

Slide the reservoir selector to the right.

Connect the Arm.

2




Ensure both of the white clips on the arm tubes are open.

Close the valve on the pressure bulb. Pressurise the fluid unit into the green zone.

IMPORTANT

Do not overpressurise (into the grey zone) as this will cause blood to leak from the puncture sites under the Skin.

3

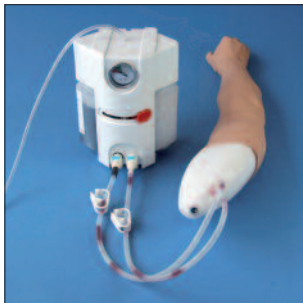

Air pressure in the right reservoir forces the blood out through the Arm and into the left reservoir.

Eventually the right reservoir will empty and air will flow through the veins. Once the amount of blood coming out is minimal, release the pressure and disconnect the Arm.

Emptying the fluid unit prior to storage/transport

1



IMPORTANT
When emptying the fluid unit the blood to be removed must always be in the left hand reservoir. See 'Emptying the Veins' (page 7).

Slide the reservoir selector to the left.

Connect the refill bottle.

2



Close the valve on the pressure bulb. Pressurise the fluid unit into the grey zone.

Air pressure forces the blood from the left reservoir into the bottle.

Keep the pressure in the grey zone to ensure that the bottle fills quickly.

3



Once most of the blood has been removed and only air is bubbling through the left tube, release the pressure and disconnect the bottle.

Removing the Skin

1



Using both hands, roll the upper part of the Skin down over onto itself as far as it will go. Do not force the Skin.

Fold the Skin down over onto itself again.

2



Gently pull the Skin down onto the wrist.

Gently pull the Skin off the wrist and onto the hand. The hand should gradually go floppy as the Skin is worked off the arm.



Once most of the Skin is off the arm, hold the Skin by the upper part...

... and remove it.

Fitting the Skin

1



To make the Skin easier to attach, lubricate the arm with the supplied water-based lubricant.

Ensure that the hand, forearm and upper arm are well lubricated.

2



Align the Skin with the underlying anatomy of the arm.

Gently work the Skin up the arm.

3





Pull the fingers onto the underlying hand anatomy.

Slide the Skin up to the top of the arm. Ensure that the ridge around the inside engages with groove at the top.

If necessary, slide the Skin around slightly so that it fits the arm anatomy properly.

Changing the Veins



1

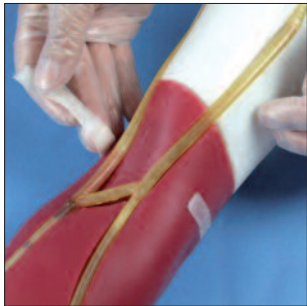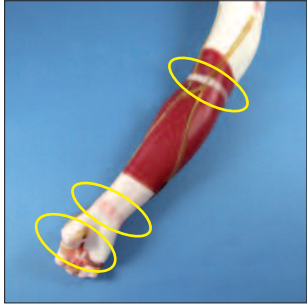

IMPORTANT
Before changing the veins follow the procedures 'Emptying the Veins' (page 7) and 'Removing the Skin' (page 9).

Remove and discard the strips of tape running around the thumb, wrist, and the upper ACF area. (New tape is supplied with the replacement veins).


2

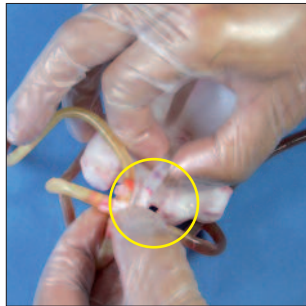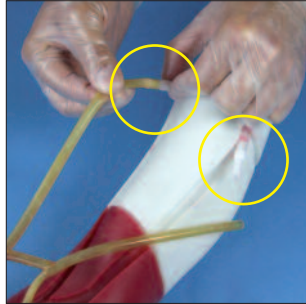

Disconnect the 2 upper connectors first.

Gently pull the lower connector out of the palm. Hold the transparent tube and pull the vein out from it.


3

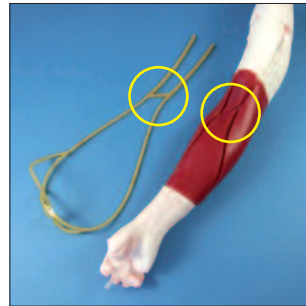

Remove the vein and discard it. Avoid spilling any residual blood left in the veins.

Lay the new vein on the work surface in the correct position (2 veins at the top, small hand veins at the bottom and the slant of the median cubital matches the groove in the vein module).


4

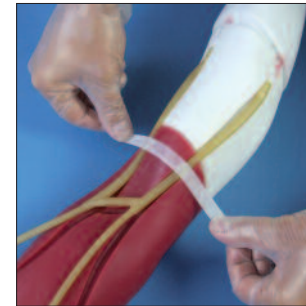

Attach the top 2 veins and push the transparent tubes back into the arm so that the connectors fit snugly.

Fit the veins into the grooves and wrap the supplied tape around the upper ACF area to hold them in place. Get another person to help you if necessary.

5

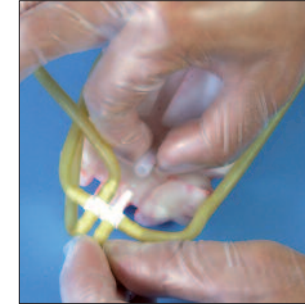

Attach the palm connector and push the transparent tube back into the arm so that the connector fits snugly.

Fit the veins into the grooves of the forearm and tape them in place around the wrist.

6



Ensure that the veins around the hand fit in the grooves correctly. Tape the vein into the thumb area.

Settle the veins into the grooves correctly. The arm is now ready to have the skin fitted.

Storing the Arm (Depending on version bought)


Cardboard box



IMPORTANT
The fluid unit should always be emptied before storing it in the cardboard box.

Follow the procedures described in 'Emptying the fluid unit' (page 8).


Case: top layer

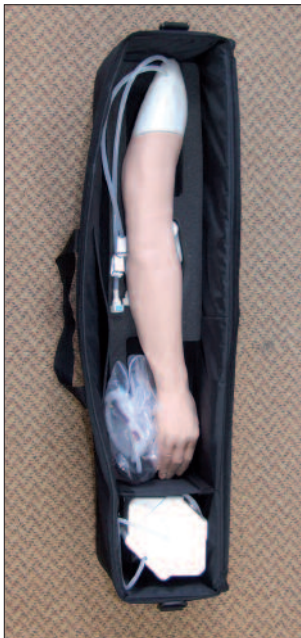

IMPORTANT
If the product is being transported long distances or being shipped, the fluid unit should always be emptied before storing it in the case.

Follow the procedures described in 'Emptying the fluid unit' (page 8).

Case: bottom layer

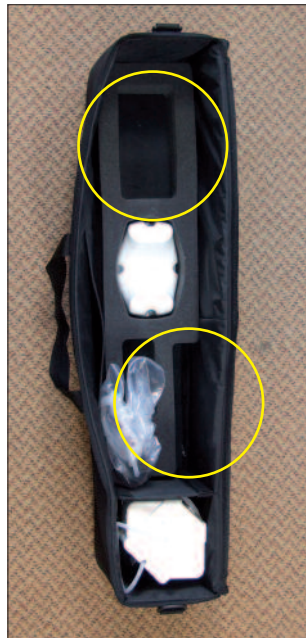

Bottom layer with additional storage for the AV Harness and additional venepuncture equipment.

Notes